

## **User Guide**

Receiver-In-Ear (RIE) Models:

Interton Centro Interton Cosmo Interton Step Interton Crisp Interton Share



Left Hearing Instrument		Right Hearing Instrument		nent
Serial number		Serial number		
Receiver type	Low Power Normal Power/Medium Power High Power Ultra Power	Receiver type	Low Power Normal Power/Medium Power High Power Ultra Power	
Reciever size	Reciever size Reciever size			
Battery size	312 brown			
Open/ standard	Small Medium Large St dome	lome Small Medium Large		Custom
fitting:	Open fitting air dome	Pow	er dome	RIE mould

Programme	Beep	Description of when to use
1	•	
2	7.7	
3	777	
4	7 7 7 7	

Note: Your hearing system might not support all 4 environmental programmes. Ask your hearing care professional for details.

### SPECIFIC FEATURES SUPPORTED BY YOUR HEARING SYSTEM:

Delayed on-activation		AutoPhone32	
Multi-Function button 24		Telecoil/Tele-loop system	
Push button	Ш	Direct Audio Input (DAI)	L
Wireless			

Ask your hearing care professional to marked options supported by your hearing system.

#### Thank You

Thank you for selecting an Interton hearing system. We are proud of our hearing products and are confident you have chosen one of the best products available.

Please familiarize yourself with the information in this guide. It contains important instructions for proper use and care, technical performance information, and other general information about your hearing system. Your hearing system has been adjusted to your particular hearing loss. Your hearing health care professional will explain these adjustments and the special features of your particular model.

### **Becoming Accustomed to Amplification**

While purchasing a hearing system is a major step, it is only one step in a process toward more comfortable hearing. Successfully adapting to the amplification your hearing system provides takes time and consistent use.

You will enjoy more benefits from your Interton hearing system by taking the following actions:

- Wear the system regularly in order to get comfortable with using it.
- It takes time to get used to a hearing aid. It may help to begin by wearing your hearing aid for short periods even as little as 15 minutes and then gradually building up your wearing time. In a way, it's no different than adjusting to contact lenses. Speak to you hearing care professional, who can design a schedule tailored just for you.
- As you get more comfortable with the system, increase the wearing time and wear your hearing system
  in multiple types of listening environments.

It may take as long as several months for your brain to get used to all the "new" sounds around you. Following these suggestions will give your brain time to learn how to interpret amplification and increase the benefits you get from using an Interton hearing system.

Hearing instrument type designations for models included in this user guide are Model: SY312, FCC ID: X26SY312, IC: 6941C-SY312, Model: MRIE, FCC ID X26MRIE, IC: 6941C-MRIE, Model SY312e, FCC ID: X26SY312e, IC: 6941C-VE312 and Model VE312, FCC ID: X26VE312, IC: 6941C-VE312. Please see pages 12 and 14 for list of models referring to all types.

#### Statement:

This device complies with Part 15 of the FCC Rules and IC rules.

Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules and ICES-003 of the IC rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates,

uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from the one in which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Changes or modifications can void the user's authority to operate the equipment.

#### Intended use

Generic air-conduction hearing instruments are wearable sound-amplifying devices intended to compensate for impaired hearing. The fundamental operating principle of hearing instruments is to receive, amplify, and transfer sound to the ear drum of a hearing impaired person.

The products are in compliance with the following regulatory requirements:

- In EU: the device conforms to the Essential Requirements according to Annex I of Council Directive 93/42/EEC for medical devices (MDD) and essential requirements and other relevant provisions of Directive 1999/5/EC (R&TTE).
- The declaration of conformity may be consulted at www.interton.com
- In US: FCC CFR 47 Part 15, subpart C.
- Other identified applicable international regulatory requirements in countries outside EU and US. Please refer to local country requirements for these areas.

- In Canada: these hearing instruments are certified under the rules of IC.
- Japanese Radio Law and Japanese Telecommunications Business Law Compliance. This device is granted pursuant to the Japanese Radio Law (電波法) and the Japanese telecommunications Business Law (電気通信事業法) This device should not be modified (otherwise the granted designation number will become invalid)
- Patents: US 7,593,537 US 8,000,849

#### Contents


Hearing instrument(s) model and basic	Push button (optional)26	Care and maintenance39, 40	Troubleshooting guide
information2, 3	Flight mode (wireless devices only)28, 29	Cleaning the reciever tubes and domes 42	Technical Data60
Thank you	Telephone use30	Cleaning RIE moulds	Warranty and repairs 64
Becoming Accustomed to Amplification 5	Listen to radio or TV	Cleaning wax guard	Temperature test, transport and
Intended use	Using Hearing Instruments with iPhone®,	How to apply domes	storage information65
Instruments component identification 13	iPad®, and iPod touch® (Interton Centro 6)31	Mini domes	Hearing instrument expectations66
On/Off function	Cellular phones	Standard domes	Warning to Hearing Aid Dispensers
Delayed activation (optional) 16	AutoPhone (optional) 32-36	Using Interton hearing instrumetns with	(US only)
Inserting/ Replacing the battery 17	Telecoil (optional)	smart phone apps	Important Notice for Prospective
Low battery indicator	Tele-loop systems (optional)	Wireless warnings and precautions	Hearing Aid Users (US only) 67
Inserting/removing hearing instruments 20	Direct Audio input (optional)	(wireless devices only)	Children with hearing loss (US only) 69
Multi function button (optional) 24	Connecting/Disconnecting audio boots 38	General warnings49	

Mico RIE hearing instruments are available in the following variants: C662-DR. Cl362-DR

Mini RIE hearing instruments of type SY312 are available in the following variants: C661-DRW, C461-DRW

Mini RIE hearing instruments of type MRIE are available in the following variants:
CI360-DVIRW, CI260-DVIRW, ST660-DRW, ST460-DRW
ST260-DRW


- 1. Receiver tube
- 2. Receiver Open Dome
- 3. Receiver Tulip Dome
- 4. Receiver Power Dome
- 5. RIE mould
- 6. Sports lock7. Battery compartment
- 8. Receiver
- 9. Push button
- 10. Multi-Function button
- 11. S receiver tube
- 12. NP receiver tube
- 13. HP receiver tube
- 14. HP2 receiver tube
- 15. Left/right indicator
- 16. Model, and serial number
- 17. Direct Audio Input

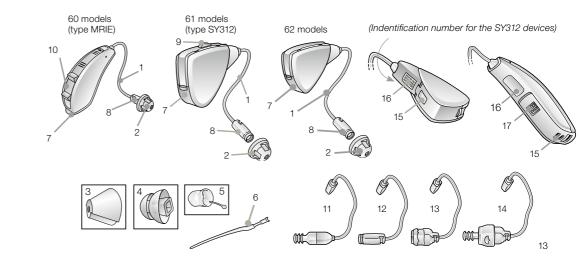

Mini Receiver In-the-Ear (RIE) hearing instruments of type VE312 with FCC ID: X26VE312, IC number 6941C-VE312 and size 312 battery are available in the following variant: CT661-DRW, CT461-DRW, and CT261-DRW.

The identification number for the VE312, instrument models can be found at location "15" as indicated in the illustrations on page 15.

- Receiver tube
- 2. Receiver Open Dome
- 3. Receiver Tulip Dome
- 4. Receiver Power Dome
- RIE mould
- 6. Sports lock
- 7. Battery compartment (serial number inside)
- 8. Receiver
- 9. Push button
- 10. LP receiver tube
- 11. MP receiver tube
- 12. HP receiver tube
- 13. UP receiver tube
- 14. Left/right indicator
- 15. Model name and serial number

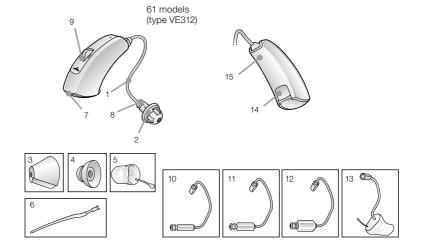


### Getting started

#### On&Off function

- 1. When the battery door is closed, the hearing instrument turns on, and the default program will be activated.
- 2. To turn off the hearing instrument, open the battery door. Many individuals can use their fingernail to pull it open.

Whenever the hearing instruments are not in use, remember to turn them off to avoid unnecessary battery consumption.

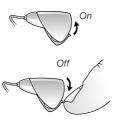

### Delayed activation

Hearing instruments can be turned on once you have placed them on your ears. If you prefer to turn them on just prior to placing them on your ear, your hearing care professional can activate a function called Delayed on-activation. This function will delay the time in which the hearing instruments turn on by several seconds after the battery compartment is closed. With Delayed on-activation, a beep will be heard for each second of the delay period.

### Inserting/Replacing the battery

- 1. Open the battery door completely by using your fingernail.
- 2. Remove the used battery if present. Insert the new battery with the positive side in the correct position.
- Gently close the battery door.



Always use new Zinc-Air batteries that have a minimum remaining shelf-life of 1 year.





- Remove the batteries to prevent leakage when the hearing instruments are not in use for an extended period of time.
- 2. Do not attempt to recharge batteries (Zinc Air) which are not specifically designated as rechargeable because they may leak or explode.
- 3. Do not place batteries in your mouth. Consult a physician immediately if a battery has been swallowed, as they can be harmful to your health.
- 4. Keep batteries away from pets, children and individuals who are mentally challenged.
- Do not attempt to dispose of batteries by burning them. Used batteries are harmful to the environment. Please dispose of them according to local regulations or return them to your hearing care practitioner.

#### Low battery indicator

Your hearing care professional can set your hearing instrument to give an acoustical indication when the battery is reaching its end of life. The hearing instrument will reduce amplification and emit a melody if battery power gets too low. This signal will recur every five minutes until the hearing instrument automatically switches off. It is recommended that you keep spare batteries on hand.

### Low battery indicator (instruments paired with accessories only)

Active usage of the Interton wireless accessories (Remote Control, Phone Clip and TV Streamer) requires more battery power from the hearing instruments than when these are working on their own. When the battery in the hearing instrument has depleted to a level at which use of the Interton TV Streamer TV and Phone Clip cannot be supported, the hearing instrument will play two sets of descending tones. After this, your hearing instrument and Interton Remote Control will continue to work as usual, but you will not be able to use your Interton TV Streamer and Phone Clip. At some point the battery level will not support the remote control either and you will once again hear the descending tones. The hearing instruments will continue to work as usual. Once a new battery is inserted, full operation of the accessories will resume.

### **Sports lock**

Sports lock will be applied or adjusted by your hearing care professional.

### Inserting/Removing hearing instruments

#### Insertion (custom RIE moulds)

- Hold the RIE mould between your thumb and index finger and position its sound outlet in your ear canal.
- 2. Slide the RIE mould all the way into your ear with a gentle, twisting movement.
- 3. Move the RIE mould up and down and gently press to ensure it is positioned correctly in the ear. Opening and closing your mouth can ease insertion.
- 4. Make sure the hearing instrument is seated behind the ear..

By experimenting, an easier method may be discovered. With proper insertion, hearing instruments should fit snugly but comfortably. If the hearing instruments cause irritation of the ears, contact your hearing care professional.

## ↑CAUTION

Never attempt to modify the shape of the hearing instrument, earmoulds, or tubing yourself.

It may be helpful to pull your ear up and outward with your opposite hand during insertion.



#### Insertion (domes)

- 1. Hold the receiver tube where it bends, and gently push the dome into the ear canal. Push the dome far enough into the ear canal so that the receiver tube lies flush with the side of the head.
- 2. It is important that the tube and the dome fit correctly into your ear.
- 3. When the dome is placed correctly, you should not be able to see the receiver tube sticking out when facing a mirror.

### **CAUTION**

You should never attempt to bend or modify the shape of the thin tube.

## **A**CAUTION

Use only original Interton/GN Hearing consumables e.g. tubes and domes.



#### Removal (RIE molds)

- 1. Grasp the removal string and pull the RIE mould outward.
- 2. Consult your hearing care professional if you have difficulties removing the hearing instrument.

#### Removal (domes)


- Hold the receiver tube with your thumb and forefinger and pull the tube outward.
- 2. Consult your hearing care professional if you have difficulties removing the hearing instrument.



### **Operation of the hearing instrument**

#### Multi-Function button (optional)

The multi-function button is designed to change the volume or listening programs of the hearing instrument, based on different ways it is pressed.

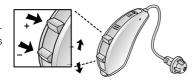


When using the program or multi-function buttons to switch programs, each press will move the instrument to the next program. For example, if it was in program 1 it will switch to program 2, if it was in program 2 it will switch to program 3 etc. When you close the battery door and switch the instrument on, it will start in program 1. Press the program or multi-function buttons if you want to move to a different listening program.

If necessary, your hearing care practitioner can change these settings and fill in the following table to indicate the new settings:


Multifunction button action	Default setting	New setting
Short press up	Increases volume	
Short press down	Decreases volume	
Long press up (2 seconds)	Changes programs	
Long press down (2 seconds)	Activates streaming	

#### Push button (optional)

Depending on your experience level with hearing instruments, individual hearing needs, and the type of listening environments you experience, your hearing care professional may activate additional programmes. If additional programmes have been activated, the following list explains how they work.

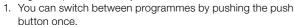

- 2. You will then hear one or more beeps. The number of beeps indicates which programme you have selected (one beep = programme one, two beeps = programme two, etc.).
- 3. When the hearing instruments are turned off and then back on, the hearing instrument always returns to the default setting (programme one).

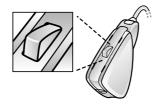

#### Only applicable to wireless devices.

If your hearing system supports wireless functionality and it is paired with wireless accessory like TV Streamer you can activate streaming mode.

- 1. Push and hold the push button or multi-function button down for 2 seconds.
- 2. You will then hear short melody that indicates streaming mode.

To switch back to environmental program push the push button shortly. Hearing Instrument will return to the default setting (programme one).

For easier everyday use of your wireless hearing instrument controls you can use wireless remote control. Ask your hearing care professional for more information.

### Flight mode. Only applicable to wireless devices\*



When boarding a flight or entering an area where RF transmitters are prohibited, wireless functionality must be deactivated, as it is not allowed to radiate radio signals during flights or in otherwise restricted areas.

#### For Crisp devices:

It is possible to disable wireless operation by opening and closing the battery compartment of the hearing instrument while at the same time pressing the push button. When disabled manually, wireless operation may be re-enabled by opening and closing the battery compartment normally, (i.e. without at the same time pressing the push button).

\*Instruments paired with accessories only

#### For Centro and Step devices:

It is possible to disable wireless operation by opening and closing the battery compartment three times within a ten second period (open-close, open-close, open-close). Your instruments will now be in flight mode.

If the hearing instrument is in flight mode, the hearing instrument must have been operating in flight mode for at least 10 seconds before attempting to enable wireless again. it is possible to re-enable wireless operation by opening and closing the battery door once. 10 seconds after this operation is completed, wireless operation will begin again.

It is important to wait an additional 15 seconds after wireless function resumes before opening and closing the battery compartment again for any reason. If the battery compartment is opened and closed during this 15 second window, flight mode will resume.

#### Telephone use

If your hearing instruments are fit with a receiver open dome or receiver tulip dome, you can probably use the telephone as you normally would by holding it up to your ear canal opening. If your hearing instruments are fit with a receiver power dome or RIE mould, finding the optimal position for holding a telephone while using a hearing instrument may require practice for some individuals, and one or more of the following suggestions may be helpful.

- 1. Hold the telephone as you would normally.
- 2. Hold the telephone towards the top of the ear (closer to where the microphones are).
- 3. If whistling occurs, it may take a few seconds of holding the telephone in the same position before the hearing instrument eliminates the feedback.
- Any whistling may also be decreased by holding the telephone slightly away from the ear.





5. Depending on your individual needs, your hearing care professional may activate a programme specifically for telephone use.

### Listen to radio or TV

When listening to the TV or the radio, start out by listening to news commentators since they usually speak clearly, then try other programmes. If you find it difficult to listen to TV or radio, your hearing care professional will be able to give you advice on available accessories to enhance your listening capabilities for TV and radio.

# Using Hearing Instruments with iPhone®, iPad®, and iPod touch® (Interton Centro 6)

Interton Centro 6 is a Made for iPhone instrument and allows for direct communication and control with an iPhone, iPad, or iPod touch. For assistance in pairing and using these products with your Interton Centro 6 device, please contact your hearing care professional.

### **Cellular phones**

Your hearing instrument is designed to comply with the most stringent Standards of International Electromagnetic Compatibility. However, not all cell phones are hearing instrument compatible. The varying degree of disturbance can be due to the nature of your particular cellular phone or of your wireless telephone service provider.

If you find it difficult to obtain a good result while using your cellular phone, your hearing care professional will be able to give you advice on available accessories to enhance listening capabilities.

### AutoPhone (optional)

The AutoPhone function, allows your hearing instrument to automatically switch to your telephone programme when a telephone receiver with an AutoPhone magnet is raised to the ear. When the telephone receiver is removed from the ear, the hearing instrument automatically returns to the previous listening programme.





#### Placement of AutoPhone magnets

Place AutoPhone magnet on your telephone receiver to allow operation of the AutoPhone function. In order to place AutoPhone magnet properly:

- 1. Clean the telephone receiver thoroughly.
- 2. Hold the telephone vertically, in a position similar to when making a telephone call.
- 3. Place the magnets just below the telephone receiver. Make sure not to cover the microphone openings. If necessary, move the magnet to another position to improve ease of use and comfort while speaking.
- 4. If you are not satisfied with the strength of AutoPhone, you can reposition the AutoPhone magnet or add additional AutoPhone magnets.

Only use a recommended cleaning agent to clean the telephone prior to placing the magnet on the phone.

#### AutoPhone usage

Telephones can be used in a normal manner. A short melody will indicate that the AutoPhone feature has automatically switched the hearing instrument to your telephone programme. Initially, you may need to move the telephone receiver slightly to find the best position for reliable AutoPhone activation and good hearing on the telephone. When remove you the telephone receiver AutoPhone will stay activated for a few seconds to avoid accidental switching off. After that hearing instrument will switch to previously used environmental programme.



#### AutoPhone warnings

- 1. Keep magnets out of reach of pets, children and mentally challenged persons. If a magnet is swallowed, please seek advice from a medical practitioner.
- The AutoPhone magnet may affect some medical devices or electronic systems. The manufacturer of any magnetically sensitive devices (e.g. pacemakers) should advise you regarding appropriate safety precautions when using your hearing instrument and magnet in close proximity to the medical device or electronic system in question.
- 3. if the manufacturer cannot issue a statement, we recommend keeping the magnet or a telephone equipped with the magnet 30 cm (12") away from magnetically sensitive devices (e.g. pacemakers).



AutoPhone precautions

- High distortion during dialing or phoning may mean that the magnet is not in the optimal position relative to the telephone receiver. To avoid the issue, please move the magnet to another place on the telephone receiver.
- 2. Only use magnets supplied by Interton/GN Hearing.

### Telecoil (optional)

If equipped, a telecoil can be activated by your hearing care professional and accessed through one of the additional programmes. A telecoil picks up a telephone's magnetic signal and converts it to sound. An optional telephone programme may help to improve speech understanding on the telephone. When using a telecoil programme, the receiver of the telephone may need to be held closer to the hearing instrument. The handset of the telephone may need to be moved to slightly different positions in order to find the best reception.

### Tele-loop systems (optional)

Many places, such as theatres, houses of worship, and schools are equipped with tele-loop systems. When using a telecoil programme with tele-loop systems, sound is picked up directly and may improve speech understanding. If there is no sound from the hearing instruments in a tele-loop system and with a telecoil programme activated, the tele-loop system may not be turned on or is not operating correctly. If a facility is not equipped with a tele-loop system, sitting as close as possible to the front may be helpful.

### **Direct Audio Input (optional)**

Use of Direct Audio Input (DAI), which enables a direct connection of the hearing instruments to items such as television, radio, and remote microphones, may increase speech understanding for some individuals. The sound source is connected to the hearing instruments by a cable or a wireless FM system to the audio boot. This accessory connects to the bottom of the hearing instruments, and once properly clicked into place, the hearing instruments switch to DAI automatically.

#### Connecting/Disconnecting audio boots

Connecting audio boots

 Align the tip of the audio boot with the groove just above the battery compartment and below the model number.



- 2. Once in place, move the boot in the direction of the battery compartment. (2)
- 3. Gently click the audio boot onto the hearing instrument.

### Disconnecting audio boots


- 4. Press and hold the button on the front side of the audio boot.
- 5. Gently remove the audio boot from the hearing instrument.

### Care and maintenance



Please follow the following instructions to prolong the durability of your hearing instruments:

 Keep your hearing instrument clean and dry. Wipe the case with a soft cloth or tissue after use to remove grease or moisture. Do not use water or solvents, as these can damage the hearing instrument(s).

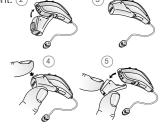

- 2. Never immerse hearing instruments in water or other liquids, as liquids may cause permanent damage to the hearing instruments.
- 3. Avoid rough handling of hearing instruments or dropping them on hard surfaces or floors.
- 4. Do not leave hearing instruments in or near direct heat or sunlight, such as in a hot, parked car, as excessive heat can cause damage or deform the casing.
- 5. Do not wear your instrument while showering, swimming, in heavy rain or in a moist atmosphere such as a steam bath or sauna.
- 6. If your instrument does get wet, or if it has been exposed to high humidity or perspiration, it should be left to dry out overnight with the battery out and the battery compartment open. It is also a good idea to put the instrument and battery in a sealed container together with a drying agent (desiccator) overnight. Do not use the instrument until it is completely dry. Consult your hearing care professional as to which drying agent to use.
- 7. Remove your hearing instrument when applying such things as cosmetics, perfume, aftershave, hair spray, and suntan lotion. These might get into the instrument and cause damage.

#### The receiver tube

The receiver tube for most Interton products comes in three different power levels:

Small (S), Normal Power (NP) and High Power (HP/HP2). Depending on your hearing loss you have either been fitted with the S, NP or the HP/HP2 receiver type.

For Interton Centro and Interton Step 6, the receiver tube comes in four different power levels: Low Power (LP), Medium Power (MP), High Power (HP) and Ultra Power (UP).

The receiver tube contains the wiring to the receiver which delivers the sound to the ear canal. It is important that the receiver tube and the receiver dome/RIE mold fits correctly in your ear. If the receiver tube or the receiver dome/RIE mould irritates your ear in any way and prevents you from wearing your hearing instrument, please contact your hearing care professional. You should never attempt to modify the shape of the receiver tube yourself. The receiver tube and the receiver dome/RIE mould should be cleaned regularly. Please see instructions in the next section.

#### Cleaning the receiver tubes and domes

The receiver tube and the receiver dome should be cleaned regularly. Use a damp cloth to clean the receiver tube and receiver dome on the outside. Do not use water when you are cleaning the receiver tubes or the receiver domes. Please see instruction on page 41 or 43 for how to change the wax guard filter.

Note: Thin tube and dome systems should be changed every three months or sooner, should the components become stiff or brittle.


### **Cleaning RIE moulds**

- 1. Separate the mould from the receiver tube.
- Clean the RIE mould using a mild soap, and rinse with lukewarm water.
- 3. After cleaning, dry RIE moulds thoroughly and remove any residual water and debris from the tubing utilising an air bulb and wire loop.




### Changing wax guard for NP receiver tube

- 1. Clean any debris from the old waxquard.
- 2. Insert the wand into the old waxquard.
- 3. Twist the wand with the waxquard in a clockwise direction to ensure it is attached to the wand.
- 4. Pull the wand and waxquard away from the tube/mould.
- 5. Insert the old waxguard into the center of the HF3 wheel.
- 6. Dispose of the old waxguard by drawing the wand to the narrow end of the center disposal area.

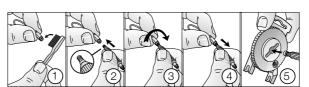

- 7. Insert the empty wand into a new waxguard on the HF3 wheel.
- 8. Pull the new waxguard attached to the wand away from the HF3 wheel.
- 9. Insert the wand into the receiver tube/mould.
- 10. Twist the wand to release the new waxguard onto the receiver tube/mould.

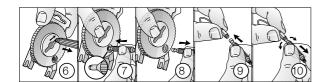


### Changing wax guard for HP receiver tube

#### Power Dome:

Please contact your hearing care professional to have the wax guard in the power dome exchanged for you.

### RIE Mould:

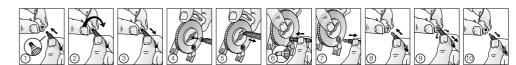

### How to apply domes

It is recommended that your hearing care professional change domes, as incorrect dome replacement could result in the dome falling out in the ear.

#### Mini domes


- 1. Push the new dome over the flanges on the receiver.
- 2. Make sure that the new dome is properly and securely mounted.



#### Standard domes

Standard domes are mounted in a similar manner to mini domes, but a few extra steps are required. Tulip domes consist of two "petals". It is important to note that the largest petal is the outermost petal. To ensure this:

- 1. Push the largest petal away from the receiver using a finger. This bends the petal forward.
- 2. Then push the largest petal back, and it will be placed on top of the smaller petal.

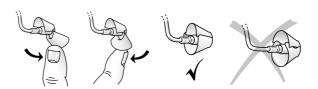

### Using Interton hearing instruments with smart phone apps

1

Intended use of smart phone apps:

Interton smart phone apps are intended to be used with Interton wireless hearing aids. Interton smart phone apps send and receive signals from the Interton wireless hearing aids via smart phones for which the apps have been developed.

#### Use with smart phone apps:

- Notifications of app updates should not be disabled, and it is recommended that the user installs all updates to ensure that the app will function correctly and will be kept up to date.
- The app must only be used with Interton devices for which it is intended, and Interton takes no responsibility if the app is used with other devices.
- If you would like a printed version of the user guide for a smart phone app please consult customer support or our website to obtain a printed user guide.

## Only applicable to wireless devices $\bigwedge$

**Depoy:** CAUTION eneral precaution (wireless hearing instrume

General precaution (wireless hearing instrument): When the wireless function is activated, the device uses a low-powered digitally coded transmissions in order to communicate with other wireless devices. Although unlikely, nearby electronic devices may be affected. In that case, move the hearing instrument away from the affected electronic device. When using wireless functionality and the devices are affected by electromagnetic interference, move away from the source of interference.

## **⚠** WARNING

General warnings (hearing instrument): Be careful when boarding flights, to remember to deactivate the wireless functionality. Turn off your wireless functionality by using the flight mode in areas where radio frequency emission is prohibited.

For use of wireless functionality only use Interton/ GN Hearing accessories. For further guidance regarding e.g. pairing, please refer to the user guide of the relevant Interton/ GN Hearing accessory.

This device operates in the frequency range of 2.4 GHz - 2.48 GHz.

This device includes a RF transmitter that operates in the range of 2.4 GHz - 2.48 GHz.

### Applicable to all devices



General warnings

- 1. Consult a hearing care professional if you discover a foreign object in your ear canal, if you experience skin irritation, or if excessive ear wax accumulates with the use of the hearing instrument.
- 2. Different types of radiation, e.g. from NMR, MRI or CT scanners, may damage the instrument. Therefore, do not wear the instrument during these or other corresponding scanning procedures. Other types of radiation (burglar alarms, room surveillance systems, radio equipment, mobile telephones etc.) will not damage the instrument. They could, however, momentarily affect the sound quality or create strange sounds from the instruments.
- 3. Do not wear hearing instruments in mines, oil fields, or other explosive areas unless those areas are certified for hearing instrument use.

- 4. Do not allow others to use your hearing instruments. This may cause damage to the hearing instruments or to the hearing of the other individual.
- Instrument usage by children or mentally challenged persons should be supervised at all times to
  ensure their safety. The hearing instrument contains small parts that could be swallowed by children.
  Please be mindful not to leave children unsupervised with this hearing instrument.
- 6. Hearing instruments should be used only as prescribed by your hearing care professional. Incorrect use may result in sudden and permanent hearing loss.
- 7. External devices connected to the electrical input must be safe according to the requirements of IEC 60601-1-1, IEC 60065, or IEC 60950-1, as appropriate.

## **♠** CAUTION

Only connect Interton hearing instruments to Interton/GN Hearing accessories intended and qualified to be used with Interton hearing instruments.

If device is broken, do not use.

SYMPTOM	CAUSE	POSSIBLE REMEDY	PAGE IN THIS USER GUIDE
No sound	Not turned on	Turn on by closing the battery door	16
	Dead battery	Replace battery	17
	Battery door will not close	Insert battery properly	17
	Blocked RIE mould or dome	Clean RIE mould or dome	42
	Blocked wax filter	Replace wax filter or consult your hearing care professional	33, 43

SYMPTOM	CAUSE	POSSIBLE REMEDY	PAGE IN THIS USER GUIDE
Not loud enough	Incorrect RIE mould placement	Reinsert RIE mould	20
	Blocked RIE mould or dome	Clean RIE mould, replace dome, replace filter	42, 43, 44
	Change in hearing sensitivity	Consult your hearing care professional	-
	Excessive ear wax	Consult your hearing care professional	-
	Volume set too low	<ul> <li>Increase the volume control if available or consult your hearing care professional</li> </ul>	24, 25

SYMPTOM	CAUSE	POSSIBLE REMEDY	PAGE IN THIS USER GUIDE
Excessive whistling / feedback	Incorrect RIE mould placement	Re-insert RIE mould carefully	20
	Incorrect dome placement	Re-insert dome	22
	Excessive ear wax	Consult your hearing care professional	-
	Feedback control may need adjustment	Consult your hearing care professional	-
	RIE mould worn or damaged	Consult your hearing care professional	-
	Hearing instrument settings not optimal	Consult your hearing care professional	-

SYMPTOM	CAUSE	POSSIBLE REMEDY	PAGE IN THIS USER GUIDE
Sound distorted / not clear	Weak battery	Replace battery	17
	Improper earmould or dome fit	<ul> <li>Consult your hearing care professional</li> </ul>	-
	Hearing instrument damaged	Consult your hearing care professional	-
	Hearing instrument settings not optimal	Consult your hearing care professional	-
Wireless does not work	Possible Root Cause - Device is in flight mode	<ul> <li>For Crisp devices with push button: Open and close the vices without push button: Open and close the battery</li> <li>For all Centro and Step devices: Open and close the bareactivate 10 seconds later. (If Root Cause is device in f</li> </ul>	door twice within 10 seconds attery compartment once. Wireless will

If there are any other problems not mentioned in this guide, please contact your hearing care professional.

Techr	Technical Data				
Small/Low receiver	Hearing Instrument Model	Maximum output (2ccCoupler / IEC 60118-7)	Battery size		
	ST460 S	114 dB SPL (typical)	312		
	ST660-DRW, CT661-DRW, CT461-DRW, CT261-DRW	113 dB SPL (typical)	312		

Tec	Technical Data				
Power	Hearing Instrument Model	Maximum output (2ccCoupler / IEC 60118-7)	Battery size		
Normal/Medium Por	C461-DRW NP Open, C661-DRW NP Open, C1360-DVIRW NP Open, C1260-DVIRW NP Open, SR1360-DVIR NP Open	114 dB SPL (typical)	312		
Norm	ST660-DRW, CT661-DRW, CT461-DRW, CT261- DRW	116 dB SPL (typical)	312		

Tec	Technical Data				
al Power receiver	Hearing Instrument Model	Maximum output (2ccCoupler / IEC 60118-7)	Battery size		
	C461-DRW NP, C661-DRW NP, CI360-DVIRW NP, CI260-DVIRW NP, SR1360-DVIR NP	114 dB SPL (typical)	312		
	ST460-DRW NP, ST260-DRW NP	114 dB SPL (typical)	312		
Normal	C662-DR NP Open, Cl362-DR NP Open	114 dB SPL (typical)	10A		
	C662-DR NP, CI362-DR NP	114 dB SPL (typical)	10A		

Ta alami	! D-+-
Lechni	cal Data

High Power Receiver	Hearing Instrument Model	Maximum output (2ccCoupler / IEC 60118-7)	Battery size
	C461-DRW HP, C661-DRW HP, CI260-DVIRW HP, CI360-DVIRW HP, SR1360-DVIR HP	117 dB SPL (typical)	312
	CI362-DR HP, C662-DR HP	117 dB SPL (typical)	10A
	ST460-DRW HP, ST260-DRW HP	119 dB SPL (typical)	312
	ST660-DRW, CT661-DRW, CT461-DRW, CT261-DRW	122 dB SPL (typical)	312
Ultra Power Receiver	ST660-DRW, CT661-DRW, CT461-DRW, CT261-DRW	129 dB SPL (typical)	312

### Warranty and repairs

Interton provides a warranty on hearing instruments in the event of defects in workmanship or materials, as described in applicable warranty documentation.

In its service policy, Interton pledges to secure functionality at least equivalent to the original hearing instrument.

As a signatory to the United Nations Global Compact initiative, Interton is committed to doing this in line with environment-friendly best practices. Hearing instruments therefore, at Interton's discretion, may be replaced by new products or products manufactured from new or serviceable used parts, or repaired using new or refurbished replacement parts.

The warranty period of hearing instruments is designated on your warranty card, which is provided by your hearing care professional.

For hearing instruments that require service, please contact your hearing care professional for assistance.

Interton hearing instruments that malfunction must be repaired by Interton qualified technician. Do not attempt to open the case of hearing instruments, as this will invalidate the warranty.

### Temperature test, transport and storage information

Interton Hearing Instruments are subjected to various tests in temperature and damp heating cycling between -25° C (-13F) and +70° C (+158F) according to internal and industry standards.

During transport or storage, the temperature should not exceed the limit values of  $-20^{\circ}$  C (-4F) to  $+60^{\circ}$  C (+140F) and relative humidity of 90% RH, non condensing (for limited time). The air pressure between 500 and 1100 hPa is appropriate.



### Hearing instrument expectations

A hearing aid will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions

Consistent use of the aid is recommended. In most cases, infrequent use does not permit you to attain full benefit from it.

The use of a hearing aid is only part of hearing rehabilitation and may need to be supplemented by auditory training and instructions in lip-reading.



### **⚠** Warning to Hearing Aid Dispensers (US only)

A hearing aid dispenser should advise a prospective hearing aid user to consult promptly with a licensed physician (preferably an ear specialist) before dispensing a hearing aid, if the hearing aid dispenser determines through inquiry, actual observation, or review of any other available information concerning the prospective user, that the prospective user has any of the following conditions:

- (i) Visible congenital or traumatic deformity of the ear.
- (ii) History of active drainage from the ear within the previous 90 days.

- (iii) History of sudden or rapidly progressive hearing loss within the previous 90 days.
- (iv) Acute or chronic dizziness.(v) Unilateral hearing loss of sudden or recent onset within the previous 90 days.(vi) Audiometric air-bone gap equal to or greater than 15 decibels at 500 hertz (Hz), 1,000 Hz. and 2.000 Hz.
- (vii) Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.
- (viii) Pain or discomfort in the ear.

### Important Notice for Prospective Hearing Aid Users (US only)

Good health practice requires that a person with a hearing loss have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before purchasing a hearing aid. Licensed physicians who specialize in diseases of the ear are often referred to as otolaryngologists, otologists or otorhinolaryngologists. The purpose of the medical evaluation is to assure that all medically treatable conditions that may affect hearing are identified and treated before the hearing aid is purchased. Following the medical evaluation, the physician will give you a written statement that states that your hearing loss has been medically evaluated and that you may be considered a candidate for a hearing aid. The

physician will refer you to an audiologist or a hearing aid dispenser, as appropriate, for a hearing aid evaluation. The audiologist or hearing aid dispenser will conduct a hearing aid evaluation to assess your ability to hear with and without a hearing aid. The hearing aid evaluation will enable the audiologist or dispenser to select and fit a hearing aid to your individual needs. If you have reservations about your ability to adapt to amplification, you should inquire about the availability of a trial-rental or purchase-option program. Many hearing aid dispensers now offer programs that permit you to wear a hearing aid for a period of time for a nominal fee after which you may decide if you want to purchase the hearing aid. Federal law restricts the sale of hearing aids to those individuals who have obtained a medical evaluation from a licensed physician. Federal law permits a fully informed adult to sign a waiver statement declining the medical evaluation for religious or personal beliefs that preclude consultation with a physician. The exercise of such a waiver is not in your best health interest and its use is strongly discouraged.

### Children with hearing loss (US only)

In addition to seeing a physician for a medical evaluation, a child with a hearing loss should be directed to an audiologist for evaluation and rehabilitation because hearing loss may cause problems in language development and the educational and social growth of a child. An audiologist is qualified by training and experience to assist in the evaluation and rehabilitation of a child with hearing loss.

Be aware of information marked with the warning symbol



**WARNING** points out a situation that could lead to serious injuries, **CAUTION** indicates a situation that could lead to minor and moderate injuries.



Advice and tips on how to handle your hearing instrument better.



Equipment includes RF transmitter.

Interton Centro 6 is compatible with iPhone 6, iPhone 6 Plus, iPhone 5s, iPhone 5c, iPhone 5, iPad Air 2, iPad Air, iPad (4th generation), iPad mini 3, iPad mini 2, iPad mini with Retina display, iPad mini and iPod touch (5th generation) using iOS 7.X or later. Apple, the Apple logo, iPhone, iPad and iPod touch are trademarks of Apple Inc., registered in the U.S. and other countries.

"Made for iPhone" means that an electronic accessory has been designed to connect specifically to iPhone and has been certified by the developer to meet Apple performance standards. Apple is not responsible for the operation of this device or its compliance with safety and regulatory standards. Please note that the use of this accessory with iPhone may affect wireless performance.







Please ask your local hearing care professional concerning disposal of your hearing instrument



Any issues relating to the EU Medical Device Directive 93/42/EEC, or Council Directive 1999/5/EC on Radio Equipment and Telecommunications terminal equipment should be directed to Interton A/S.

Worldwide headquarters Interton A/S | Lautrupbjerg 7 | DK-2750 Ballerup | Denmark Tel.: +45 45 75 1111 | Fax: +45 45 75 1119 | www.interton.com | CVR no. 55082715