



BTE MODEL

## For a life on the GO

Intuitive and wireless, Interton Step hearing aids bring back the sounds your customers enjoy, in every situation. Innovative features such as Personal Noise Reduction, Environmental Gain Tuner and Ear-to-ear communication available in a wide variety of models make Interton Step perfectly suited for up to 90% of people with hearing loss. With compatibility to the new 2nd generation Wireless Accessories you can offer a complete hearing system that gives your customers an enjoyable life.


## Standard Configuration

- Connectivity to wireless accessories
- Wireless ear-to-ear communication
- Nano Coating technology
- Supports Open and Closed configurations
- Push button and programmable Volume Control
- Size 13 battery
- Battery door with integrated On/Off switch
- 5 different colors
- Telecoil with T and MT modes
- Direct Audio Input
- Supports various domes and ear hooks for comfortable fit

## Fitting Requirements

- Interton Appraise 2.4 Fitting Software (or higher)
- Programming adaptor with CS44 socket cable
- Fitting interface (see table)

## Fitting Range

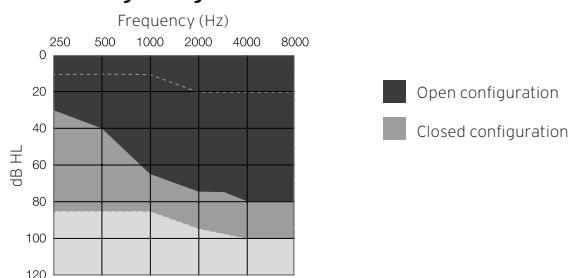

## Model Overview

ST670-DW  
ST470-DW  
ST270-DW

SOUND PROCESSING	Step 6	Step 4	Step 2
Logarithmic WDRC, number of channels	17	14	12
Gain Handles in FSW	9	7	6
Total comfort programs	4	4	3
<b>BASICS</b>			
Signal tones	●	●	●
Power-on Delay	●	●	●
<b>NOISE REDUCTION</b>			
Microphone Noise Reduction	●	⊙	○
Adaptive Noise Reduction	●	⊙	○
Personal Noise Reduction	●	●	-
Adaptive Wind Noise Reduction	●	⊙	○
<b>SPEECH UNDERSTANDING</b>			
Integrated Directionality	●	-	-
Fixed Directionality	●	●	●
Omni Directionality	●	●	●
Speech-Focused Directionality	●	●	●
Auto-Steered Directionality	●	●	●
Selectable Beam Width	●	●	○
SoundScape technology	●	●	●
<b>FEEDBACK MANAGEMENT</b>			
Feedback Manager Plus	●	●	●
Feedback Manager Plus - Music Mode	●	●	-
Pre-set Feedback Manager	●	●	●
Fitting Protection	●	●	●
<b>PROTECTION</b>			
Nano coating	●	●	●
<b>SPECIALS</b>			
Datalogging	●	●	●
AutoPhone	●	●	●
Environmental Gain Tuner	●	⊙	○
In-situ	●	●	●
<b>2.4 GHZ WIRELESS TECHNOLOGY</b>			
ear-to-ear communication	●	●	●
Volume Control synchronization	●	●	●
Push Button synchronization	●	●	●
Comfort Phone	●	●	●
Wireless Phone Clip	●	●	-
Wireless Remote Control	●	●	●
Wireless TV Streamer	●	●	-
<b>FITTING INTERFACE</b>			
Wireless fitting with Airlink	●	●	●
Hi-Pro	●	●	●
NOAHlink	●	●	●
SpeedLink	●	●	●

○ Standard

⊙ Advanced

● Ultimate

## Technical specifications

		Open		Closed		
		IEC 60118-0 IEC 711 Ear simulator	IEC 60118-7 ANSI S3.22 2cc coupler	IEC 60118-0 IEC 711 Ear simulator	IEC 60118-7 ANSI S3.22 2cc coupler	
Reference test gain (60 dB SPL input)	1600 Hz/HFA	44	41	47	47	dB
Full-on gain (50 dB SPL input)	Max.	62	51	66	59	dB
	1600 Hz/HFA	53	46	56	52	
Maximum output (90 dB SPL input)	Max.	130	122	134	128	dB SPL
	1600 Hz/HFA	126	118	130	123	
Total harmonic distortion	500Hz	0.2	0.1	0.5	0.4	%
	800 Hz	0.5	0.2	0.6	0.4	
	1600 Hz	0.4	0.5	0.6	0.5	
Telecoil sensitivity (1 mA/m input)	Max.	90	-	96	-	dB SPL
Full-on Telecoil sensitivity @ 1mA/m	1600 Hz / HFA	82	75	85	81	
HFA - SPLIV @ 31.6 mA/m (ANSI)	HFA	-	100	-	107	
Equivalent input noise w/o noise reduction		23	23	23	22	dB SPL
1/3 Octave equivalent input noise, w/o noise reduction	1600 Hz	10	-	10	-	
Frequency range (DIN 45605)		100-6900	100-6770	100-7080	100-6850	Hz
Current drain (quiescent / operating)		1.1 / 1.2	1.1 / 1.2	1.1 / 1.1	1.1 / 1.2	mA

Data in accordance with IEC 60118-0, IEC 60118-7, Supply Voltage 1.3 V.

Patents pending

All specifications are subject to change without notice

400287000-GB-16.02-Rev.B

